# INFRARED

## BODY THERMOMETER

**NON-CONTACT**

PITHG02

USER MANUAL

platinet

## EN User manual

### Non-contact electronic infrared body thermometer

This manual is part of the product and should be stored in a way that allows you to read it in the future. It contains important information about the settings and operation of the device.

The set includes:

- Electronic thermometer equipped with an infrared sensor
- User manual

Precautions:

Before using the thermometer for the first time, please read the following notes:

- This product is only intended for the use described in this manual
- Measurement performed with this product cannot substitute for a proper medical examination
- The infrared sensor is the most delicate part of the thermometer. Do not touch the lens of the infrared sensor with your fingers.
- If the product is not used for more than two months, remove the batteries and store them separately.
- Dispose of used batteries in accordance with applicable law. Batteries may explode if overheated. After discharged, they must be collected separately from other waste.
- Do not submerge the product in water or expose it to sunlight, which may damage it.
- Do not drop or expose to electrical surges.
- Do not use outdoors or near a fireplace or stove.
- If necessary, clean the surface of the glass with the infrared sensor with 70% isopropyl alcohol.
- The device should be used in an ambient temperature between 5 ° C and 40 ° C (0 ° F to 140 ° F). The thermometer should not be used outside this range. To get the most accurate measurement results, the thermometer and the test person should be in the same room for about 15-20 minutes to equilibrate their temperature before the measurement.
- Maintenance tasks that you can perform only include battery replacement and cleaning.

Thermometer application:

This product measures body temperature by measuring the infrared radiation emitted from the skin surface on the forehead. The non-contact thermometer is designed to measure body temperature (Body mode) and surface temperature (Surface mode). It can also check the temperature of the room, food, bath water, baby milk, etc.


Product description and method of operation:

The thermometer consists of a plastic housing, LCD display, function buttons, infrared temperature sensor, components and electronic circuits.

The non-contact thermometer is a handheld, battery-powered device designed to periodically measure and monitor the temperature of the human body by measuring the infrared radiation emitted from the skin surface on the forehead of infants, children and adults. The optical system of the thermometer directs and focuses the infrared energy to the detector, and then the electronics of the instrument convert the signal from the detector to the reading shown on the display.


Drawing and description of the device structure:

1. Infrared sensor
2. Battery socket
3. LCD display
4. "+" button
5. "-" button
6. "M" button, ON / OFF
7. Temperature measurement button

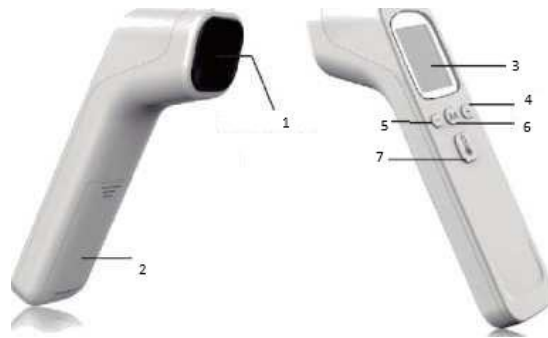Display drawing and description:

1. Surface mode - measuring surface temperature
2. Body mode - measuring body temperature
3. Measurement result
4. Number of the saved measurement
5. Icon of speaker and volume level
6. Measurement result in ° C scale
7. Measurement result in ° F scale
8. Low battery indicator

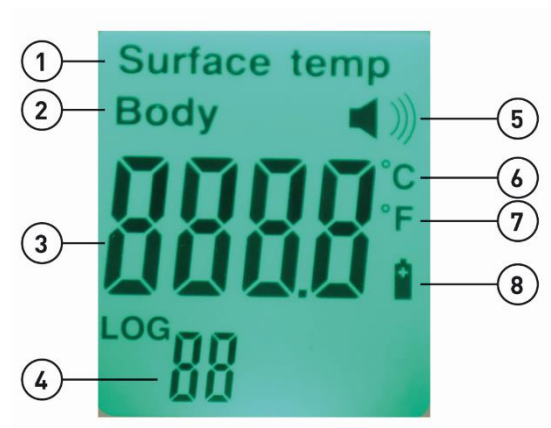Taking measurements

Comments:

- Measurements should be made only on the exposed part of the forehead, any covering (i.e. hair or a hat) may cause an erroneous measurement.
- Accurate measurement requires cleaning the test surface (i.e. forehead) from make-up and other substances that could distort the measurement result. Bring any hair off your forehead, wipe away any sweat.
- If you are testing the temperature of babies who are usually wrapped in a blanket or clothes, the neck should be exposed for at least 1 minute before measuring.

- Cold compresses and other methods of reducing fever may affect the correct course of the measurement and its results.
- Measurements should be made in an environment with a stable temperature distribution. Do not measure near an operating fan or air conditioner, in a strong blast of cold air.
- Do not take measurements in the area of strong sunlight.
- Hold the thermometer approximately 3-5 cm from your forehead, press the ON / OFF scan button. The temperature will be displayed immediately.
- If the measurement result seems to you inadequate to the actual body temperature you feel, repeat the measurement with the sensor pointing at the skin behind the ear.

The measuring procedure:

How to turn on the device?

- Press the "M" button to start the device
- Press the "M" button again to turn off the device.

The thermometer automatically turns off after 60 seconds of non-use

How to measure the temperature?

- Turn on the device
- Aim the meter towards your forehead at a distance of 3-5 cm.
- Press the largest button with a thermometer symbol, to measure the temperature.

How to switch the temperature measurement mode?

- Turn on the device.
- Press the "-" button to switch the body temperature measurement mode to surface or object temperature measurement.

How do I change the temperature unit of measure from ° C to a ° F?

- Turn on the device.
- Press and hold for 5 sec. "-" button, until "F1" appears in the left corner of the display.
- Press the button 3 times and change the symbol "F1" to "F4".

The symbol appears in the right corner of the display temperature measurement units.

- Change between ° C / ° F by pressing the "+" button.
- Turn off the device to save changes.

How to enter the history of temperature measurements?

- Turn on the device.
- Press the "+" button to view the history of results temperature measurements.
- Press the "+" button again to view other saved data.
- Turn off the device.

How do I silence the audible signal?

- Turn on the device.
- Press and hold the "-" button for 5 seconds, until "F1" appears in the left corner of the display.
- Press the "+" button to turn the sound on / off (the volume symbol means that the sound is on)
- Turn off the device to save changes.

Attention:

If the measurement result seems to you inadequate to the actual body temperature you feel, repeat the measurement with the sensor pointing at the skin behind the ear.

Cleaning and storage:

- After use, wipe the thermal sensor with a damp cotton swab, dry with a paper towel and store in the package.
- The device should be cleaned once a week with a cloth soaked in water or a neutral detergent, avoiding to flood it or excessively wet it.
- The device should not be repaired by yourself. Do not modify this thermometer without the approval of the manufacturer. Do not disassemble or attempt to repair the device or its component yourself.

## Battery change

The device requires two AAA (LR03) batteries. Open the battery compartment. To do this, slide down the battery compartment cover on the back. Take out the used batteries from the battery compartment. Install new batteries. Make sure the batteries are inserted the right way round. Close the battery compartment.

## Security

Temperature measurement results should not be used as the basis for medical examination. For measuring babies, keep the baby still. Movement can cause erroneous readings. It is not allowed to modify or disassemble the device. If the thermometer is not used for an extended period of time, remove the batteries and store in a safe place. The device is not suitable for use in the presence of flammable anesthetics or oxygen. This device has been tested and found to comply with the limits for a medical device in EN 60601-1: 2006 + A11: 2011 + A1: 2013 + A12: 2014. The limits are designed to provide reasonable protection against harmful interference in a typical medical installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to other equipment in the vicinity. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause interference to other equipment, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Change the direction or position of the receiving device.
- Increase the separation between the devices.
- Connect the equipment into an outlet on a circuit different from that to which the other equipment or devices are connected.
- Consult the manufacturer's service technician for assistance. Do not measure during or immediately after breastfeeding. Temperature readings may vary by location. Therefore, we recommend that you measure the temperature at 2 points. Do not use the thermometer in environments with high humidity.

## Technical parameters:

Product name: non-contact infrared thermometer

Operating ambient temperature: 16 °C -35 °C

Working atmospheric pressure: 70kPa-106kPa

Working air humidity: ≤85%

Supply voltage: 3.0 V DC (2 AAA alkaline batteries)

Size: 140 \* 40 \* 40mm

Weight: 71 g / 88.5 g

Measuring range: Body temperature: 32 °C - 42.9 °C

Measurement accuracy: 35.0-42.0 °C:  $\pm 0.2$  °C; outside the range of 35.0-42.0 °C:  $\pm 0.3$  °C

Measurement distance: 3-5 cm

Measurement time: <2 seconds

3-color display: body temperature: 34.0-37.3 °C green; 37.4-38 °C yellow; 38.1-42.6 red

Automatic shutdown: 15 seconds

Storage and transportation: -20 °C -55 °C, Relative Humidity:  $\leq 93\%$

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