CORDLESS MASSAGE GUN

PMAS1

USER MANUAL

 platinet


Round Head

Used for hand, back, leg, muscle.

Bullet Head

Used for deep tissue, such as joint, palm and plantar.

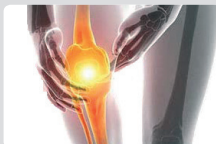Flat Head

Used for all parts of the body.

„U” Shape Head

Used for neck, spine, and Achilles.

EN User Manual

Package Contents
Massage Gun
4pcs Interchangeable Heads
USB Charging Cable
Instruction Manual

Key Features

Rechargeable Cordless Massage
6 Power Level Massage Mode
Helps relieve muscle soreness and stiffness
Improves range of motion
Promotes blood circulation

Safety Precautions

For adults only.
Do not place anything into the ventilation holes on the shaft.
Do not put this product into water or let water get into it.
Use only the attached USB Charging Cable for recharging.
Do not leave this product unattended when working or on charge.
Do not use accessories which are not recommended.
Do not use this product outdoors.
Recommended usage time of this product is 15 minutes/session.
Please do not use it in a moist environment such as bathroom.
Please do not use it when there is a drastic change in temperature.
Please switch off the power and unplug the machine after use.
Please do not use it when surrounding temperature is above 104°F (or 40°C).
Please do not use it when surrounding temperature is below 41°F (5°C).
Do not place it near fire or directly under the sun.

People who have the below symptoms, please consult a doctor before use

People who are suffering from skin conditions such as swollen skins, sunburnt or other dermatology diseases.
People who are suffering from malignant tumors, acute diseases, heart diseases, stroke or severe hypertension.
People who are suffering from any blood circulation disorders such as thrombosis, severe aneurysms and venous tumors; or from any skin infections.
People who are suffering from osteoporosis.
People who are suffering from heart disease or carrying electronic medical equipment such as heart pacemaker.
Pregnant women or women that are in menses period.
People whose bodies are wet are forbidden to use this product.
People who have open wounds.

Safety

Use only on a dry surface.
Check the massaged area frequently and immediately

stop using when there is pain or discomfort.
This product is not to be used for medical treatment purposes.

Recommended usage time of this product is 15 minutes/session.

Do not use wet hands to pull the plug.
If the power cord is damaged, it must be replaced.
If the product is found damaged, broken, suffering from power leakage, or with exposed parts, please stop using it.

Do not turn on the product when it is being charged.
Do not use this product continuously for over 30 minutes, or else the temperature control will automatically turn off the machine. If so, please allow a 30-minute break before resuming.

How to Operate

Choose a suitable massage head and plug into the massage gun.
Press the Power [Mode] button to turn on.
Switch to 6 different massage mode. *Auto Off Setting at around 15 minutes
Set the massage at a desired speed on your body, apply pressure accordingly.

How to Charge

Before using it for the first time, charge it fully (about 3 hours).
Connect the USB Charging Cable (included) to a suitable USB power outlet rated DC 5V/2A, and then the other end to the DC Charging Port of the device.
6 Red LEDs will flash, when the battery is being charged.
All 6 LED Lights will stay on, when the battery is fully recharged.

Maintenance and cleaning

Please don't disassemble the product.
Please always turn off the main power after use.
Keep device in dry and dustless environment.
Do not keep this product at a high temperature or near fire. Do not keep under direct sunlight.
Please clean this product with dry cloth.
Before cleaning, please make sure the switch has been turned off and the power plug has been disconnected from the product.

Product Specifications

Product Weight: 1237 grams
Charging Time: About 3 Hours
Working Time: About 30Minutes
Speed at Mode 1: 2100rpm
Speed at Mode 2: 2300rpm
Speed at Mode 3: 2500rpm
Speed at Mode 4: 2700rpm
Speed at Mode 5: 2900rpm
Speed at Mode 6: 3100rpm
Static State Current [Max.]: 30uA
Working Current: [Max.]: 4000mA
Static State Power [Max]: 0.22mW
Voltage: 7.4V 1500mAh
Rated Power: 20W
Rated time: 15 minutes

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PLATINET S.A.